

General Meditech 2800
G3 Intelligent General Monitoring System







General Meditech 2800




G3 Intelligent General Monitoring System

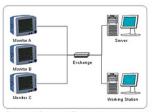

Connected System by Ethernet line, the signal is stable and reliable

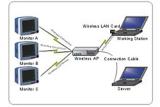

Establish LAN with AP, all monitors of the system are connected by wireless way to achieve the flexible effect.

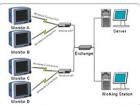

Adopt multi-point connection mode to enlarge the cover range of wireless signal. Especially when there is more sickrooms or the sickrooms are located in different floor. It's a good solution to the matter of multi floors and large quantity monitors to set distributed signal area by setting wire in small area and through exchanger, then add wireless connection point according to the condition

CHARACTER

Display

- O Simultaneous display with multi 17" monitors.
- C Each screen can display 16 sets bedside monitors in three modes. Each system can connect up to 256 set bedside monitors.
- O Display of multi waveforms (ECG, RESP, PLETH, CO2, IBP etc.)and multi parameters(HR, ST, Sp02, PR, TEMP, SBP, DBP, MAP, INCO2, ETCO2).
- O More than 10,000 hours data storage is printable according to requirement. Capture and store ECG waveform.

Control

ECG lead on each monitor can be changed in central monitoring system. NIBP measurement on each monitor can be controlled by central monitoring system.

Alarm and setup

Operating menu with multi-language selection


Alarm limit, switch and status of each monitor can be setup by central monitoring system

Central monitoring system can setup its own audio and visual alarm limits.


Support long distance monitoring and diagnosis as well as software upgrade. Adaptive network connection for RS-232 or Ethernet wireless and I or wire

Three connection modes: Wire, Wireless or Multi-Point Connection Mode.



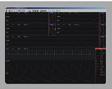



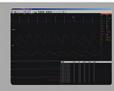

General Meditech, Inc.

Add: South office 4\F, Kezhi 1st Rd. west, Science park , Nanshan, Shenzhen, Guangdong , P.R.China Postcode: 518057

Tel: 86-755-26500832 26615545, 26546289 http://www.szmedtech.com

Fax: 86-755-26546285 E-mail:sales@szmedtech.com